



Suma Bac D10 Mini Guide - COVID-19: cleaning of non-healthcare settings



Volume of solution	500 ml bottle	750 ml bottle	1L bottle	5L container
Volume of product	10 ml	15 ml	20 ml	100 ml
Volume of water	480 ml	700 ml	1000 ml	5000 ml




Cleans and disinfects in one step



Effective against enveloped viruses according to EN14476


Use for cleaning and disinfecting:

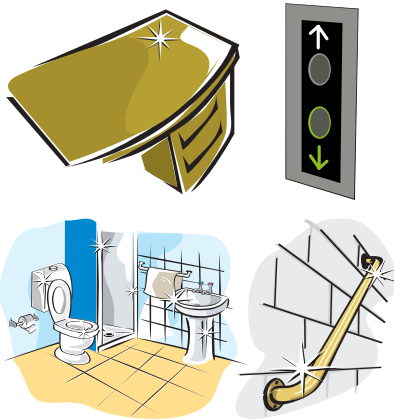

General areas: Office desks, lifts, hand rails, door handles etc.

Application...


1

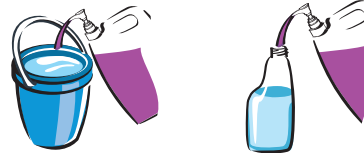

A made up solution of product will be active for up to 7 days - throw away any made up solution after 7 days. Rinse and dry bottle before refilling.


If cleaning down large, general areas it is better to use a bucket solution.

2

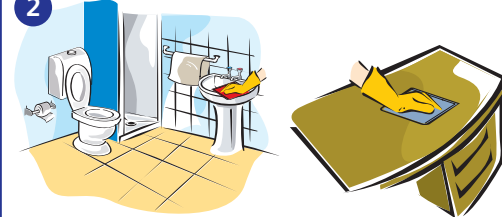

Remove debris before applying product.

3

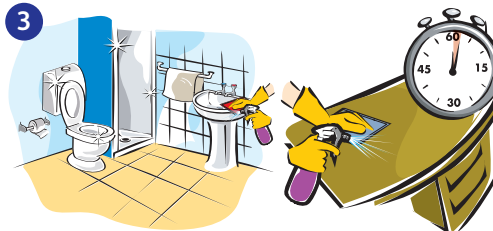

Surface must be visibly wet and allow 1 minute contact time.

4



Ensure surfaces are visibly wet after wiping. For sensitive surfaces such as lift buttons / touch pads, spray made up solution onto a cloth and wipe, then re-wipe.

www.gov.uk/government/publications/covid-19-decontamination-in-non-healthcare-settings

Ensure all surfaces are rinsed at the end of the day

